
# **Cordless Water Flosser Kabellose Munddusche**

**Model | Modell: EWF-32**

**User Manual  
Bedienungsanleitung**

**EN** **DE**

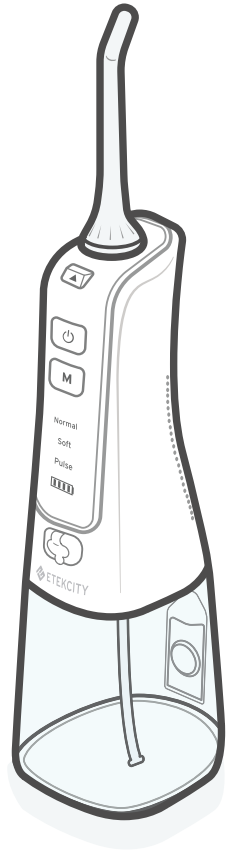**Questions or Concerns? • Fragen oder Probleme?**  
[support.eu@etekcity.com](mailto:support.eu@etekcity.com)

Package Contents	2
Specifications	4
Important Safeguards	6
Function Diagram	8
Controls	9
Using Your Water Flosser	9
Maintenance	10
• Charging the Battery	10
• Cleaning the Water Flosser	11
• Cleaning the Water Tank, Water Hose, and Water Filter	11
• Replacing the Jet Tip	12
• Storage	12
Troubleshooting	13
Warranty Information	14
Customer Support	14

1 x Cordless Water Flosser  
1 x Water Tank (Pre-Installed)  
4 x Standard Jet Tips  
1 x Charging Cable  
1 x User Manual

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**Français**

**Italiano**

<b>Input</b>	USB 5V, 1A
<b>Battery Type</b>	1400mAh Li-ion Battery
<b>Charging Time</b>	4 hours
<b>Water Tank Capacity</b>	300 mL / 10.1 oz
<b>Water Pulse Rate</b>	1400-1800 times/min
<b>Water Pressure</b>	276-724 kPa / 40-105 psi
<b>Dimensions</b>	7.3 x 8.2 x 30.2 cm / 2.9 x 3.2 x 11.9 in
<b>Weight</b>	0.3 kg / 0.7 lb
<b>Automatic Shutoff</b>	2 minutes

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## ENGLISH • Cordless Water Flosser EWF-32

---

**Thank you** for purchasing the **Cordless Water Flosser** by Etekcity.

If you have any questions or concerns, please reach out to our helpful Customer Support Team at **support.eu@etekcity.com**.

We hope you enjoy your new cordless water flosser!

# IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, these basic safety precautions should **always** be followed:

## READ ALL INSTRUCTIONS BEFORE USING

**DANGER:** To reduce the risk of electrocution while charging:

- **Do not** handle charging cable with wet hands.
- **Do not** submerge or drop the charging cable or water flosser in water or other liquid while it is charging.
- **Do not** use while bathing.
- **Do not** place or store the water flosser where it can fall or be pulled into a tub, shower stall, or sink.
- **Do not** reach for the water flosser if it has fallen into water. Unplug the charging cable immediately if it is plugged in.
- Check the charging cable for damage before using. **Never** use a damaged charging cable to charge the water flosser.

**WARNING:** To reduce the risk of burns, electrocution, fire, or injury to persons:


- **Do not** plug the charging cable into a voltage system that is different from the voltage system specified in **Specifications** (see page 2).
- **Do not** use this water flosser if it has a damaged cord or plug, if it is not working properly. Contact **Customer Support** (see page 14).
- **Do not** direct water under the tongue, into the ear, nose, or other delicate areas. This water flosser is capable of producing pressures that may cause serious damage in these areas.
- Use the water flosser **only** as indicated in these instructions and as recommended by your dental professional.
- Depending on your gums' health, your gums may bleed the first few times you use the water flosser. This is normal for any type of flossing. If bleeding is excessive or continues after 2 weeks, see a dental professional.
- Fill the water tank with water or other solutions recommended by dental professionals **only**.

- **Do not** drop or insert any foreign object into any opening on the water flosser.
- Keep the charging cable and water flosser away from heated surfaces.
- **Do not** use where oxygen or aerosol sprays are being used.
- **Do not** use mouthwash or rinse that contains iodine, bleach, or essential oils with the water flosser.
- **Do not** use the water flosser while wearing any oral jewellery. Remove them prior to use.
- **Do not** use if you have an open wound on your tongue or in your mouth.
- If your physician advised you to use antibiotic premedication before dental procedures, you should consult your dentist before using your water flosser or any other oral hygiene aid.
- Close supervision is necessary when the water flosser is used by children, invalids, or disabled persons.
- **Do not** allow children to play with the water flosser.
- **Only** use the provided charging cable to charge the battery.
- This water flosser contains batteries that are non-replaceable.
- Not for commercial use. Household use **only**.

**SAVE THESE  
INSTRUCTIONS**

## Function Diagram

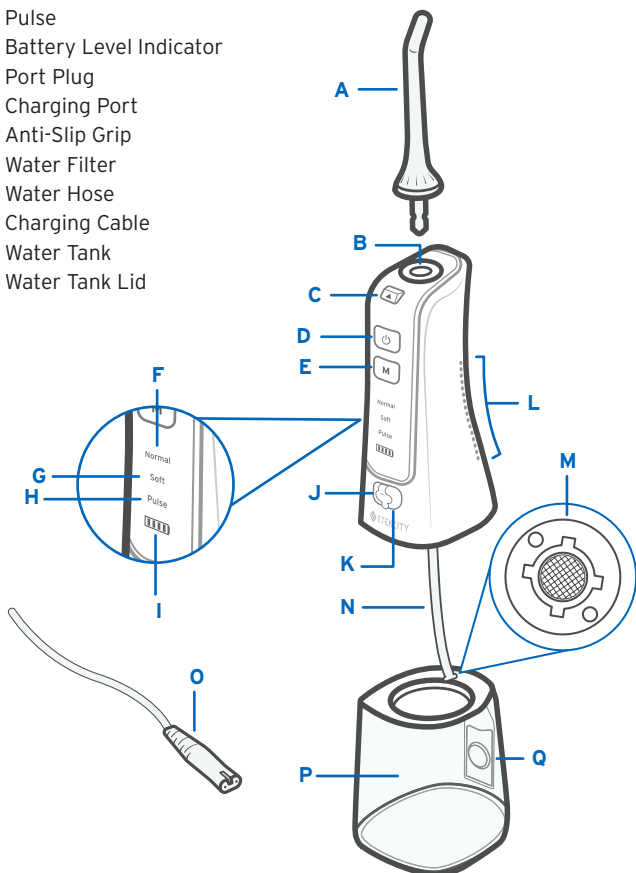- A. Jet Tip
- B. Jet Tip Port
- C. Tip Release Button
- D. On/Off Button
- E. Mode Selection Button
- F. Normal
- G. Soft
- H. Pulse
- I. Battery Level Indicator
- J. Port Plug
- K. Charging Port
- L. Anti-Slip Grip
- M. Water Filter
- N. Water Hose
- O. Charging Cable
- P. Water Tank
- Q. Water Tank Lid

## Controls

-  : Releases jet tip from water flosser.
-  : Turns the water flosser on or off.
- M** : Cycles through Normal, Soft, and Pulse modes. You can change the mode at any time, including when the water flosser is off.

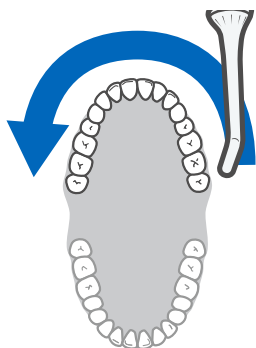
Figure 1.1

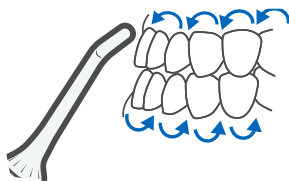
Figure 1.2

## Using Your Water Flosser

### Note:

- Before first use, fully charge the water flosser (see **Charging the Battery**, page 10).
- Use Soft mode for the first 2 weeks to allow your gums to get used to the water flosser. If bleeding continues after 2 weeks, see a dental professional.

1. Insert the jet tip securely into the jet tip port. You will hear it click into place.
2. Make sure the water flosser is off and the port plug is covering the charging port. Open the water tank lid and fill the water tank with tap water. Press the lid closed firmly.
3. Press **M** to select the water pressure mode.
4. Point the jet tip at the gumline of your back teeth and slightly close your mouth to prevent water from spraying out. Lean over a sink and press and hold  for 1 second to turn the water flosser on [Figure 1.1]. Let water flow into the sink.
5. Slide the water stream along the outside of your gums. Pause between the teeth to remove any plaque or food particles [Figure 1.2].

## Using Your Water Flosser (Cont.)

- Repeat along the inside of your gums, until all teeth are fully cleaned [Figure 1.3].
- Press  to turn the water flosser off.

### Note:

- The jet tip can be rotated 360° to adjust the direction of the water flow.
- You can switch between the 3 water pressure modes at any time.
- The water flosser will automatically shut off after 2 minutes.
- The water flosser will remember your last used water pressure.
- Water may sometimes leak from the sides of the water tank. This is normal.

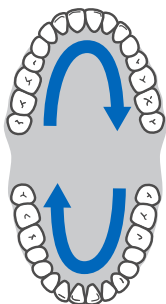
Figure 1.3

## Maintenance

### Charging the Battery

The water flosser comes with a rechargeable Li-ion battery inside. Use the included charging cable to charge the water flosser.  will turn red when it's time to recharge the battery.

- Insert the USB charging cable into the charging port on the water flosser [Figure 2.1].

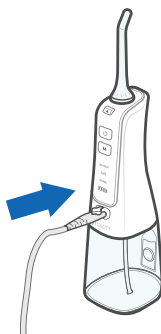
Figure 2.1

- Plug the USB charging cable into a DC 5V adapter and plug into an outlet. You can also plug the cable into a powered USB outlet.
- Charging may take up to 4 hours. The water flosser is fully charged when  turns green.

**Note:**

- **Do not** use the water flosser while it is charging.
- To maximise battery life, recharge the battery whenever  turns red.
-  indicates the water flosser has a full battery,  indicates the battery is at half power, and  indicates the water flosser must be recharged.
- The water flosser includes a lithium battery that cannot be replaced. Towards the end of the battery life, please recycle at your local recycling centre.

## Cleaning the Water Flosser

1. Turn the water flosser off.
2. Fill the water tank with clean water and allow it to run until the water tank is empty.
3. If necessary, rinse the water flosser under clean running water.
4. Dry with a soft cloth.

Figure 3.1

## Cleaning the Water Tank, Water Hose, and Water Filter

1. Turn the water flosser off and turn the water tank to the  position [Figure 3.1].
2. Gently pull off the water hose and filter from the water flosser [Figure 3.2].
3. Rinse under clean water or clean with a toothbrush or soft cloth.
4. Optionally, you can remove the water filter by gently pulling it off with a fingernail or the flat end of a toothpick [Figure 3.3]. Clean the filter with a toothbrush or soft cloth and rinse.
5. Dry all parts with a clean cloth.

Figure 3.2

Figure 3.3

## Replacing the Jet Tip

We recommend replacing the jet tip every 3 months. Press ▲ and discard the old jet tip before replacing with a new jet tip.

**Note:** To maintain the performance of your water flosser, **only** use official Etekcity jet tips. To buy replacement jet tips, go to Amazon and search [ASIN].

## Storage

1. Turn the water flosser off and empty the water tank. Wipe dry.
2. Press ▲ and pull the jet tip away from the water flosser.
3. Twist the water tank so the line is on .
4. Place the jet tip inside the water tank and screw the water tank back onto the water flosser so the line is on .
5. Store in a cool, dry area.

### Note:

- Fully charge the water flosser before storing for long periods of time.
- **Do not** store near chemicals or corrosive gases. **Do not** leave the water flosser exposed to any chemical solvent, lint, or dust.
- Keep away from heat sources and direct sunlight.

## Troubleshooting

Problem	Possible Solution
Area between the jet tip and the jet tip port leaks water.	Make sure jet tip is properly installed.
	Contact <b>Customer Support</b> (see page 14).
Water tank leaks water.	Make sure the water tank is securely in place.
Water pressure is too low.	Recharge the battery (see <b>Charging the Battery</b> , page 10).
	Clean the water flosser and jet tip (see <b>Cleaning the Water Flosser</b> , page 9).
	If the pressure is on Soft, increase the water pressure to Normal (see <b>Using Your Water Flosser</b> , page 9).
	Replace the jet tip (see <b>Replacing the Jet Tip</b> , page 12).
Water flosser doesn't turn on.	Recharge the battery (see <b>Charging the Battery</b> , page 10).
Water flosser doesn't release water.	Clean the water filter (see <b>Cleaning the Water Tank, Water Hose, and Water Filter</b> , page 11).
Water leaks from the sides of the water tank.	Because the water tank and main body of the water flosser are not completely sealed, water may leak from the sides. This is normal.

**If your problem is not listed**, please contact **Customer Support** (see page 14).

## Warranty Information

### TERMS & POLICY

Etekcitec Corporation warrants all products to be of the highest quality in material, craftsmanship, and service for 2 years, effective from the date of purchase to the end of the warranty period. Warranty lengths may vary between product categories.

Should you encounter any issues or have any questions regarding your new product, feel free to contact our helpful Customer Support Team. Your satisfaction is our goal!

## Customer Support

If you have any questions or concerns about your new product, please contact our helpful Customer Support Team.

### Customer Support

**Etekcitec Corporation**  
1202 N. Miller St., Suite A  
Anaheim, CA 92806  
USA

**Email:** [support.eu@etekcity.com](mailto:support.eu@etekcity.com)

*\*Please have your invoice and order ID ready before contacting Customer Support.*

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**Etekciti Corporation**