
## Personal Fitness Scale

Model No. EF984H

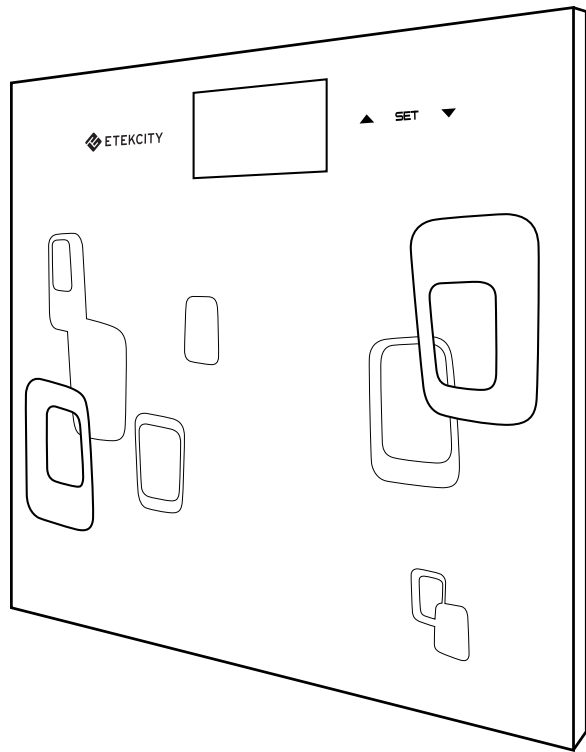EN | DE | FR

**Questions or Concerns?**  
support.eu@etekcity.com

## Thank You.

Thank you for purchasing the Personal Fitness Scale by Etekcity. This device can help you reach your physical goals by conveniently allowing you to measure body weight with reliable accuracy. Calculate body fat, water, muscle and bone percentages with the use of the BIA principal. The scale also measures suggested daily calorie intake based on personal user results. Create up to 12 individual user profiles, personalized for accuracy and reliability for the whole family. The decorative platform accents any type of home decor with a subtle modern touch. The Personal Fitness Scale is engineered to champion the concept of "smart living" through convenience.

We're constantly working at building better living with the products and services we create. Should you ever experience any issues with your new product, feel free to contact our helpful customer service team via e-mail at [support.eu@etekcity.com](mailto:support.eu@etekcity.com).

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# Table of Contents

## EN

Safety Information	6
Main Functions	7
Set-Up	8
Operation	8
Specifications	11
Customer Support	11
Troubleshooting	11
Performance Index	12
Fitness Assessment Chart	12
Warranty Information	13

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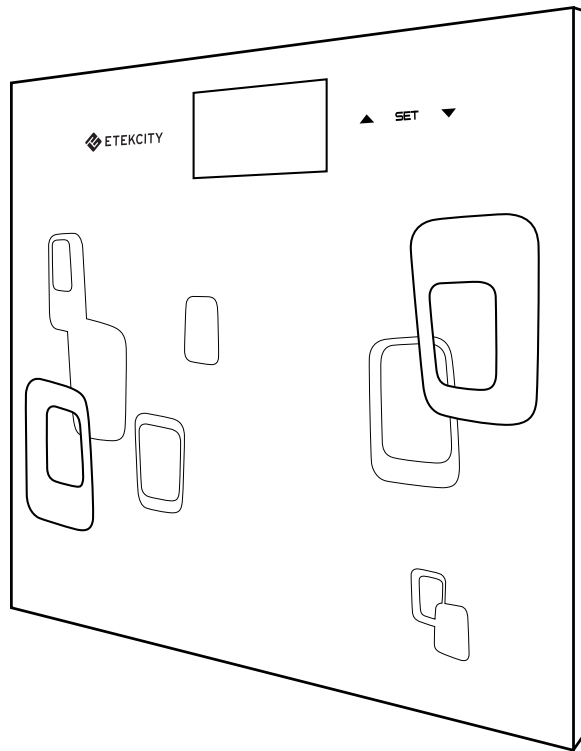
**Questions or Concerns?**  
support.eu@etekcity.com



## Safety Information

**IMPORTANT:** Please read, understand, and comply with all of the instructions provided in this manual before using the device. Failure to comply with the instructions given in this manual and/or using the device in ways other than the ones mentioned in this manual may result in serious injury and/or damage to the product.

- The glass platform is slippery when wet. Make sure the platform is dry before stepping on for measurement.
- Consult your primary-care doctor or physician before making changes to your diet, exercise plans and/or physical activities.
- This scale is not intended for commercial use.
- To weigh safely and accurately, please stand still on the center of the platform until a result is calculated.
- **DO NOT** strike or drop the scale.
- **DO NOT** stand on the edge of the scale or jump on it.
- Always place the scale on a hard, dry and flat surface before measurement.
- Handle the scale with care, as it is a precision instrument.
- Clean the scale with a damp cloth and avoid getting water or any other liquid inside the scale.
- **DO NOT** use chemical or abrasive cleaning agents to clean the scale.
- Keep the scale in a cool environment with dry air.
- Always keep the scale in a horizontal position, do not store vertically.
- The fat and hydration content evaluated by the scale are only for reference. This scale is not intended for medical purposes. If results indicate that the fat and hydration content is over or below a normal level, consult your primary physician.
- Keep the scale away from strong electromagnetic fields at all times.
- If the scale fails to turn on, check whether the batteries are installed correctly. Replace the batteries if necessary.
- If an error appears on the display or if the scale does not turn off automatically, remove the batteries for 3 seconds, then reinstall them.
- **DO NOT** overload the scale (max weight 180 kg / 400 lbs).
- If the scale is not in use for an extended period of time, remove the batteries to ensure battery conservation and to extend the scale service lifetime.
- Before measurement, always remove shoes and socks for best weight result accuracy.
- Individuals with the following conditions should not use this scale: pregnant women, those with symptoms of edema, those on dialysis treatment, and those who use heart pacemakers or other implanted medical devices.

# Main Functions

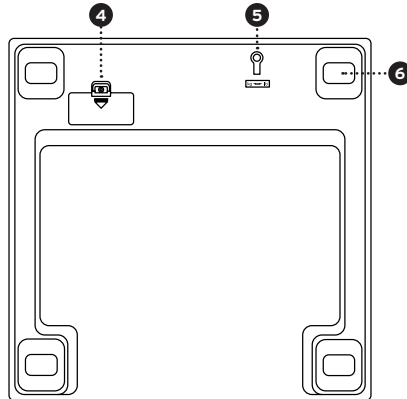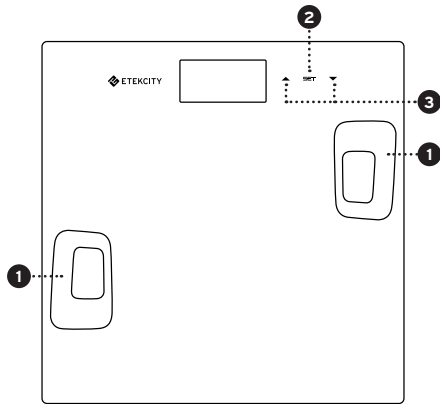
## Function Diagram

1. Electrode plates
2. Profile selection
3. Setting navigation
4. Battery compartment
5. Unit selection padding
6. Anti-skid padding

## Button Functions

- ▲ Increase the value by one interval. Press and hold to speed up value adjustment.
- ▼ Decrease the value by one interval. Press and hold to speed up the value adjustment.
- SET** Enter profiles selection and confirm selection.

## LCD Display

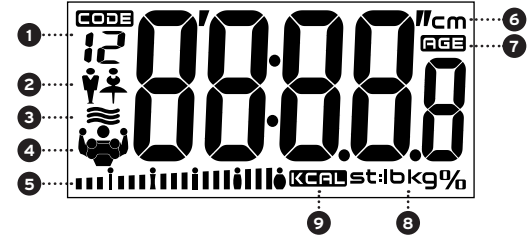
1. User memory code
2. Gender
3. Hydration %
4. Fat % | Muscle % | Bone %
5. Fat status indication
6. Height unit
7. Age
8. Weight unit
9. Calorie unit


## Set-Up

The 2 AAA 1.5V batteries are preinstalled in the scale battery compartment. Open the battery compartment cover located on the backside of the scale and remove the plastic film from the batteries.

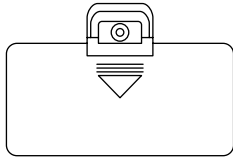
Re-apply the batteries under the correct polarity and re-seal the battery compartment in order to activate the scale.

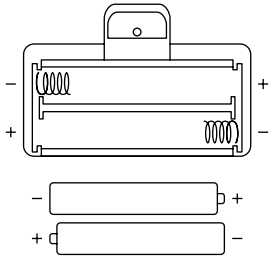
## Operation

### Normal Weighing Mode

**For best accuracy, place the scale on a hard, flat floor. Avoid carpet or soft surfaces.** Gently step on the scale, making sure your bare feet are on the **metal electrodes** on the platform surface. The scale will automatically activate.

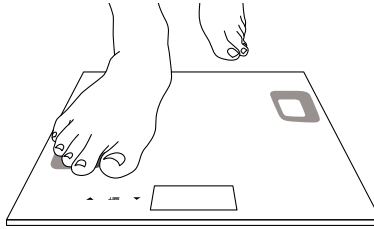Remain standing evenly on the scale, without moving or shifting your feet or body while the scale processes the results. The scale will display and lock the stabilized result.
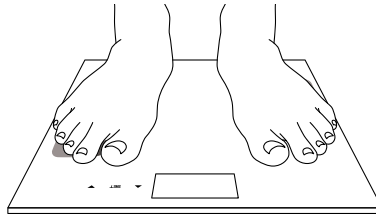
## Personal User Profiles

User Number	Gender	Height	Age
1 - 12	Male / Female	100 - 250 cm 3'0.3.5" - 8'02.5"	10 - 100

### Creating a New User Profile

1. While the scale is off, press the **'SET'** button to activate profile selection.
2. Press **'▲'** and **'▼'** to select the desired user profile number and press **'SET'** to confirm selection.
3. Press **'▲'** and **'▼'** to select the correct gender and press **'SET'** to confirm selection.
4. Press **'▲'** and **'▼'** to set your current height and press **'SET'** to confirm selection.


5. Press **'▲'** and **'▼'** to set your current age, then press firmly on the center of the platform to complete the profile customization.

◀ **NOTE:** To lock the stabilized weight results, step off the scale while the final results begin to flash.

6. The LCD screen will now display '0.0' with the profile number and gender, indicating that it is ready to calculate measurements in fitness mode.

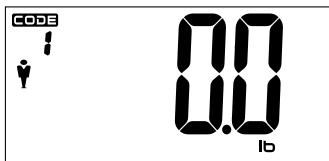
**NOTE:**

- Profile settings can be adjusted at any time by selecting the user profile and changing any value in order to maintain the most accurate results.
- If the scale is left idle for 6 seconds while in user profile settings, the changed settings will automatically save and the scale will be ready for fitness measuring mode. To return to profile settings, press and hold the 'SET' button for 3 seconds.

## Fitness Mode

**NOTE:**

- To quickly access fitness mode after weighing yourself, press and hold the 'SET' button for 3 seconds.
- Make sure to be as accurate as possible when inputting your personal profile metrics. Inaccurate profile metrics may result to inaccurate fitness readings.

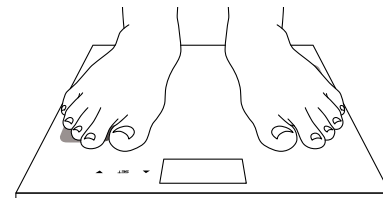1. While the scale is off, press the 'SET' button to activate profile selection.
2. Use '▲' and '▼' to select the desired user profile.
3. Once the desired profile is selected, press firmly on the center of the platform to confirm selection.
4. '0.0' will appear along with the user number and gender icon, indicating that the scale is ready for measurement use.

5. Gently step on the scale, making sure your bare feet are on the **metal electrodes** on the platform surface.
6. Remain standing evenly on the scale, without moving or shifting your feet or body while the scale processes the results. The scale will display and lock the stabilized results.

7. After displaying the calculated weight, the scale will display the results in the following order: weight, fat status and percentage, hydration percentage, muscle percentage, bone percentage, and calorie requirement. The scale will cycle through the results three times before shutting off.

Weight

Fat Status & Percentage


### Hydration Percentage

### Muscle Percentage

### Bone Percentage
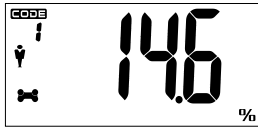
### Calorie Requirement

**NOTE:** The user calculations will be saved for reference for the next use.


## Recall Result History

1. While the scale is off, press the **'SET'** button to activate profile selection.
2. Use **'▲'** and **'▼'** to select the desired user profile.
3. Wait for 6 seconds or firmly press on the center of the platform. '0.0' will appear along with the user number and gender icon.
4. Press **'▲'** and **'▼'** to enter recall mode. The scale will cycle through the previous results 2 times before automatically shutting off.

## Converting Weight Units


Press the unit conversion button located on the backside of the scale to convert between pounds and kilograms. The LCD display will show the new unit of measurement as you switch.

**NOTE:** When setting lbs as the weight unit, the default height unit is feet/inches. When setting kg as the weight unit, the default height unit is cm.
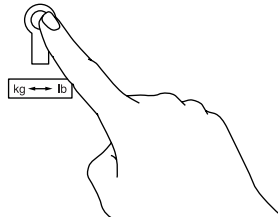
## Battery Replacement

1. Open the battery compartment cover on the backside of the scale.
2. Take out the used batteries and dispose them properly.
3. Install 2 new AAA 1.5V batteries, taking note of the correct polarity displayed in the battery compartment.
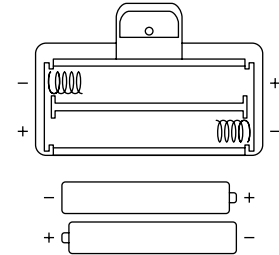
## Specifications

Fitness Metrics:	Body weight   Body fat Hydration   Muscle % Bone %
Weight Capacity:	180 kg (400 lbs)
Graduation:	0.1kg (0.2lbs)
Display:	74 x 35.6 mm
Battery:	2 x 1.5V AAA batteries
Weight Units:	kg   lb
Height Units:	cm   ft   in
User Profiles:	12
Platform:	6mm tempered glass
Dimensions:	302 x 302 x 22 mm

## Package Contents

- 1 x Personal Fitness Scale
- 1 x User Manual
- 2 x 1.5V AAA batteries (pre-installed)

## Customer Support

Should you encounter any issues or have any questions regarding your new product, feel free to contact our helpful Customer Support Team. Your satisfaction is ours!

[support.eu@etekcity.com](mailto:support.eu@etekcity.com)

## Troubleshooting

Message	Reason	Solution
"Err"	Scale overload. Maximum weight capacity: 180 kg, 400 lbs.	Step off the scale to avoid damage.
"ErrL"	Low Fat % indication.	Consult primary physician.
"Errh"	High Fat % indication.	Consult primary physician.
"Lo"	Battery power is low.	Install new batteries.
"C"	Weight reading error.	Step off the scale and retest weight.

## Performance Index

Item / Index	Body Weight	Fat %	Hydration %	Muscle %	Bone %	Calories
Index Range		4.0% ~ 60.0%	27.5% ~ 66.0%	20.0% ~ 56.0%	2.0% ~ 20.0%	
Division	0.1kg (0.2lb)	0.1%	0.1%	0.1%	0.1%	1 KCAL

## Fitness Assessment Chart

Under or at the age of 30.

Fat (F)	Hydration (F)	Fat (M)	Hydration (M)	Status	Key Line Indication
4.0 - 16.0%	66.0 - 57.8%	4.0 - 11.0%	66.0 - 61.2%	Underweight	...i
16.1 - 20.5%	57.7 - 54.7%	11.1 - 15.5%	61.1 - 58.1%	Lean	...iii
20.6 - 25.0%	54.6 - 51.6%	15.6 - 20.0%	58.0 - 55.0%	Normal	...iiii
25.1 - 30.5%	51.5 - 47.8%	20.1 - 24.5%	54.9 - 51.9%	Overweight	...iiiii
30.6 - 60.0%	47.7 - 27.5%	24.6 - 60.0%	51.8 - 27.5%	Obese	...iiiiii

Over the age of 30.

Fat (F)	Hydration (F)	Fat (M)	Hydration (M)	Status	Key Line Indication
4.0 - 20.0%	66.0 - 55.0%	4.0 - 15.0%	66.0 - 54.8%	Underweight	...i
20.1 - 25.0%	54.9 - 51.6%	15.1 - 19.5%	58.3 - 55.3%	Lean	...iii
25.1 - 30.0%	51.5 - 48.1%	19.6 - 24.0%	55.2 - 52.3%	Normal	...iiii
30.1 - 35.0%	48.0 - 44.7%	24.1 - 28.5%	52.2 - 49.2%	Overweight	...iiiii
35.1 - 60.0%	44.6 - 27.5%	28.6 - 60.0%	49.1 - 27.5%	Obese	...iiiiii

**NOTE:** Muscle mass and bone mass are not included in the assessment chart because health levels vary between individuals. For a target goal in muscle and bone mass, consult your primary health care physician.

# Warranty Information

## Terms & Policy

Etekciti warrants all products to be of the highest quality in material, craftsmanship and service for a minimum of 2 years, effective from the date of purchase. Warranty lengths may vary between product categories. For specific warranty terms, please refer to your product listing page or product operation manual.

Under the limited warranty, Etekciti will repair, replace, or refund any parts found to be defective due to manufacturer flaws. This warranty extends only to personal use and does not extend to any product that has been used for commercial use, rental use, or any other use in which the product is not intended for. There are no warranties other than that expressly set forth with each product.

This warranty is non-transferrable. Etekciti is not responsible in any way for any damages, losses or inconveniences caused by equipment failure or by user negligence, abuse, or use noncompliant with the user manual or any additional safety, use, or warnings included in the product packaging and manual.

Etekciti and its subsidiaries assume no liability for damage caused by the use of the product other than for its intended use or as instructed above and within the user manual. Some states do not allow this exclusion or limitation of incidental or consequential losses so the foregoing disclaimer

may not apply to you. This warranty gives you specific legal rights and you may also have other rights which may vary from state to state.

ALL EXPRESS AND IMPLIED WARRANTIES, INCLUDING THE WARRANTY OF MERCHANTABILITY, ARE LIMITED TO THE PERIOD OF THE LIMITED WARRANTY.

If your ETEKCITY product should require service: Should your product prove defective within the specified warranty period, please return the defective unit in its original packaging with (1) an original copy of invoice, (2) an original copy of order confirmation, (3) and your warranty ID number.

## Exceptions & Exclusions

- Damage due to abuse, accident, alteration, misuse, tampering or vandalism.
- Improper or inadequate maintenance.
- Damage in return transit.
- Unsupervised use by children under 18 years of age.

To receive a valid warranty ID number, sign into your Etekciti customer account at [www.etekcity.com](http://www.etekcity.com). At the top of the warranty page, type in your order ID number located in your order invoice to retrieve your warranty ID number. Your warranty ID number will then be sent to your email for your records. If deemed defective, your product will be repaired, replaced or refunded at no charge to you.

Date of Purchase:

---

Place of Purchase:

---

Warranty ID:

---

Order ID:

---


To register your warranty or to make a claim, please contact our helpful customer support team via email at [support.eu@etekcity.com](mailto:support.eu@etekcity.com)

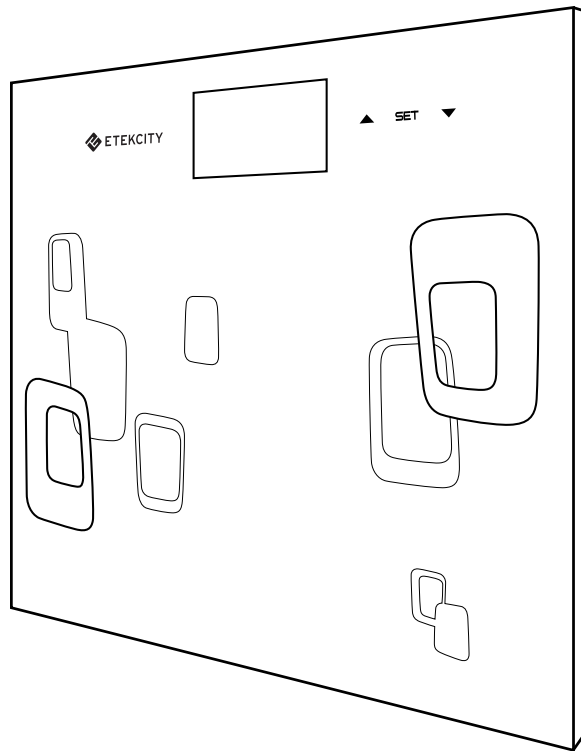*Please keep this page for your records.*

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