Shiatsu Massager

Model No.: EM-SH7


EN

DE

FR

ES

IT

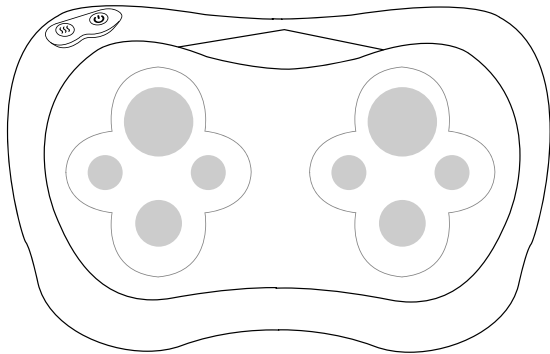EN

Function Diagram	3
Safety Information	4
Set Up	4
Operation	5
Maintenance	5
Specifications	6
Warranty	6

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



Model No.: EM-SH7

Function Diagram

1. Massager housing
2. Massager nodes
3. Power/Reverse button
4. Speed/Heat button
5. Adjustable strap
6. Power cord connector
7. AC adapter
8. Pillow Cover
9. DC car adapter

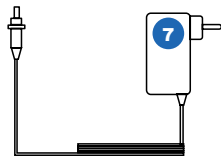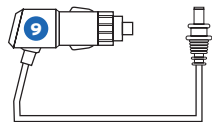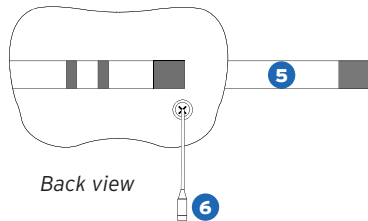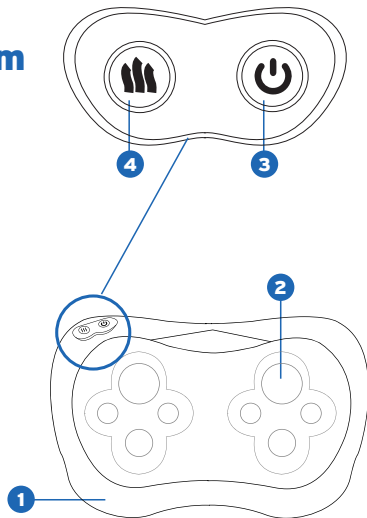Safety Information

IMPORTANT: Please read and comply with all of the instructions and warnings provided in this manual before using the product. Failure to comply with the instructions and warnings provided herein may result in inaccurate results and/or damage to the product itself.

- **DO NOT** immerse the massager in water or any other liquids.
- **DO NOT** use the device in wet or moist conditions.
- **DO NOT** operate the device if the power cord or AC adapter is damaged.
- Keep power cord and AC adapter away from heated surfaces.
- **DO NOT** use the massager under a blanket or pillow; excessive heating may occur and can possibly cause fire, electric shock, or injury.
- **DO NOT** use the massager in bed or while sleeping.
- **DO NOT** use outdoors.
- **DO NOT** stand or sit on the massager.
- **DO NOT** use metal pins or fasteners with this appliance.
- Always unplug the massager when not in use.
- Adult supervision is necessary when the device is being used near or by children.
- **DO NOT** use this product directly on swollen or inflamed areas.
- **DO NOT** use this product for more than 20 minutes at a time.
- Extensive use of the warming function may overheat the massager and shorten its life span. If the massager

overheats, discontinue use and allow it to completely cool before resuming use.

- To reduce the risk of electric shock, **DO NOT** remove the housing cover.
- Any individual who is pregnant, has a pacemaker, suffers from diabetes, phlebitis, and/or thrombosis, or is at an increased risk of developing blood clots due to recent surgery, should consult his/her physician before using a massaging device designed for home use.
- **DO NOT** use on an infant or invalid or on a sleeping or unconscious person. Do not use on sensitive skin or on a person with poor blood circulation.

Set Up

1. Remove the massager and the AC adapter from their packaging.
2. Place the pillow cover over the massager.
3. Attach the massager to an the desired chair using the adjustable strap; there are three adjustment levels for the strap.
4. Position the massager to your desired level on the chair.
5. Connect the AC adapter to the power cord connector.

NOTE: This device can be used in automobiles. Connect the DC car adapter to the power cord connector, then plug the DC car adapter into the power outlet in your car.


Operation

1. Press the Power/Reverse button to turn on the massager with heating. The massage nodes will change direction every 60 seconds.
2. Press the Power/Reverse button to reverse the direction of the massage nodes.
3. Press the Speed/Heat button once to increase or decrease the speed of the massage nodes.
4. Press and hold the Speed/Heat button for 3 seconds to turn on or turn off the warming function in the massager.
5. Press and hold the Power/Reverse button for 3 seconds to turn off the massager.

NOTE:

- The device shuts off automatically after 20 minutes.
- It is recommended to leave the device off for at least 30 minutes before the next use.

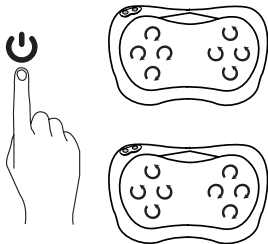Maintenance

1. Unplug the massager, allow the unit to cool before cleaning.
2. Remove the pillow cover, wash and dry it.
3. Wipe down the massager with a dry clean cloth.
4. **DO NOT** immerse the massager in water or liquid detergents.
5. If necessary, wipe down the massager with a clean, damp cloth and dry immediately with a soft, dry cloth.
6. **DO NOT** attempt to repair the massager. Contact customer service if the massager malfunctions or is in need of repair.
7. Store the massager in a cool, dry place away from sharp edges or pointed objects when not in use.

NOTE: To reduce the risk of electric shock, **DO NOT** remove the housing cover.

Specifications

AC Adapter Input	AC 100-240V / 50-60Hz
DC Car Adapter Output	DC 12V
Rated Power	24W
Size	13.3 x 4.3 x 8.6 in (340 x 110 x 220 mm)

Package Contents

- 1 x Shiatsu massager
- 1 x AC adapter
- 1 x Pillow cover
- 1 x DC car adapter
- 1 x User manual

Warranty

Terms & Policy

Etekcitec warrants all products to be of the highest quality in material, craftsmanship and service for a minimum of 2 years, effective from the date of purchase. Warranty lengths may vary between product categories.

Should you encounter any issues or have any questions regarding your new product, feel free to contact our helpful Customer Support Team. Your satisfaction is ours!

Email: support.eu@etekcity.com

*Please have your order ID number ready before contacting customer support.

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